

LipoFLUDDY staining for flow cytometer

1. Seed 300 000 cells in a 35 mm plate (or 6-well plate)
2. After 24 h add 10 μ M LipoFLUDDY. (1:1000 dilution, 1 μ L in 1mL, 1 μ L enhancer, add LipoFLUDDY and enhancer in a 1:1 ratio to eppendorf tube first, then add media and mix well).
3. Incubate the cells with LipoFLUDDY for 4hr at 37 °C in a CO₂ incubator. After incubation, wash the cells twice with 1 $\times$ PBS to remove excess dye.
4. Replace the medium containing LipoFLUDDY with normal medium or starvation medium depending on the experiment. For positive control use Torin-1 (1 - 3 μ M, maximum 12hr treatment), for negative control use Bafilomycin A1 (10 nM - 1 μ M, typically 100 nM, maximum 4hr treatment).
5. Detach the cells using 500 μ L trypsin and neutralize the trypsin with 1 mL of FACS buffer (PBS 1x + 10% serum). Spin down the cells, discard the supernatant and resuspend in 500 μ L of FACS buffer.
6. Analyse the samples using the flow cytometer.



Note

To minimize background fluorescence originating from the phenol red in the medium, it is recommended to use DMEM without phenol red or FluoroBrite DMEM.

Flow cytometer settings:

FITC = Ex: 488 nm; Em: 500-560 nm

APC = Ex: 640 nm; Em: 663-677 nm

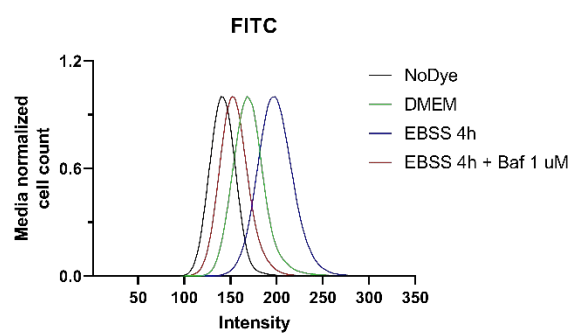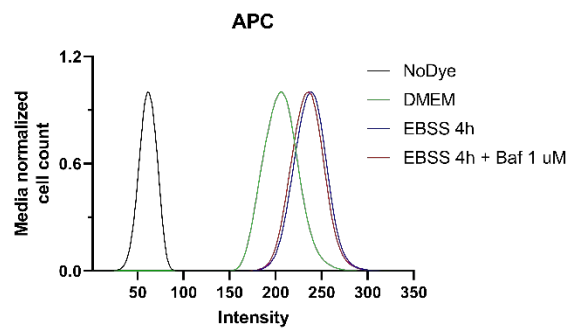

LipoFLUDDY staining for microscopy

1. Seed 300 000 cells in a 35 mm Ibidi μ -Dish.
2. After 24 h add 10 μ M LipoFLUDDY. (1:1000 dilution, 1uL in 1mL, 1uL enhancer, add LipoFLUDDY and enhancer in a 1:1 ratio to eppendorf tube first, then add media and mix well).
3. Incubate the cells with LipoFLUDDY for 4hr at 37 °C in a CO₂ incubator. After incubation, wash the cells twice with 1 $\times$ PBS to remove excess dye.
4. Replace the medium containing LipoFLUDDY with fresh medium.
5. Analyse the samples using a microscope with appropriate filters.

Note

To minimize background fluorescence originating from the phenol red in the medium, it is recommended to use DMEM without phenol red or FluoroBrite DMEM.

Microscopy settings:

Green = Ex: 488 nm; Em: 500-550 nm

Deep red = Ex: 640 nm; Em: 660-750 nm

